



wissner-  
bosserhoff  
contempora

Furniture for an ambience worth living in





“When the eyes lose their sharpness,  
the powers of understanding begin  
to increase their focus.”

*Plato*

By the year 2030 every third person in Central Europe will be over 60 years of age – a development that can already be observed today. As life span increases, society is seeing a change in values, one that offers new and exciting challenges. In order to do justice to the culture of the elderly we have to recognise and accept challenges. This is the only way how we can successfully raise the quality of living of older and especially care-needy people up to a higher level. Let's then concentrate on what is essential and do our part to help create a new care culture under a new paradigm. For we will only achieve our goals if we make humanity the focal point of both how we interact with elderly and infirm people and how we design care and nursing facilities for the elderly.



“ Each generation  
smirks about its fathers,  
laughs about its grandfathers  
and admires its  
great-grandfathers.”

*Somerset Maugham*

For residents requiring care, the third phase of life holds many changes, ones that often are accompanied by a new type of consciousness. This can find expression in one's personal environment, in a more intense way of experiencing friendship or an increased appreciation of a pleasant environment that offers individual comfort. wissner-bosserhoff is pursuing a clearly defined objective: to support the wishes and needs of people under regular care with innovative and functional products. In cooperation with sieger design, the company has now developed contempora H10, the home-like care bed for increased quality of life and comfort. We take great care to be able to provide furniture of the highest functionality that meets the demands of everyday living. At the same time, these products embody the high quality and sophisticated design for which sieger design is known. An homage to all the elderly individuals who have been willing to share their experiences and their style with us.





contemporanea H 10







01



02



03

**Previous page** contempora H10 care bed and bedside cabinet SD 2, cherry.

**01** The head board H enables the use of a homely plaid in padded look.

**02** The divided side protections lock together firmly for gapless security.

**03** The handrail provide a homely touch and safety on all sides of the bed.

**Right side** The side protections, extendable in four positions by up to 40 cm, provide individualised safety and maximum freedom of motion.



“It takes a long time before you become young.”

*Pablo Picasso*

Its bold lines and matching aesthetics allow the contempora H10 care bed to fit in harmoniously with a classic living environment. Its high-quality detail is particularly noticeable, in, for example, appealingly-shaped head and foot boards with rounded bars. The matching bedside cabinet SD 2 in cherry complements the care bed, rounding it out into a home-like unit providing individual comfort to meet the occupant's needs.

# contempora H 10





contemporanea L 10







01



02



**Previous page** Care bed contempora L10 with low foot board and bedside cabinet SD 3, maple.

**01** The divided side protections give the resident support when getting out of bed, facilitating unaided mobility.

**02** The low foot board provides an additional homely comfort and greater freedom of movement. In combination with automated comfort sit-up mechanism it allows exiting the bed from the foot end.

**Right side** Innovative elevation technology that is hidden at first look. With slender lifting column covers and the easyclick braking system.



“Everyone wants to grow old,  
but nobody wants to be old.”

*Arthur Schopenhauer*

contempora L10 complements the light, fresh atmosphere of well-lit room with its straight lines and lack of flourish. The matching bedside cabinet SD 3 in maple lends the environment a rich and harmonious look. contempora care beds are functional and yet design-orientated. Thus the innovative elevation technology, in keeping with style considerations, is modestly concealed within a slender lifting column cover. contempora L 10 is also loved for the many practical advantages it offers such as easier getting into/out of bed with its 36 cm minimum height setting. The 40 cm height adjustment range makes working with the bed especially easy on the back.

# contempora L 10



contemporanea M 09





**Previous page** contempora M09 care bed and Vivo E2 bedside cabinet, beech 01.

**01** Be on the safe side: continuous side rails railings activated by convenient pistol grip controller.


**02** Indirect lighting – The night light function aids residents in getting their bearings, even when the hour is late.

**03** Convenient: The bedside cabinet SD3 with integrated extending bedside table.

**Right side** The bedside cabinet SD2 without bedside table. In combination with a server, this variation is a solution that is both economic and individually usable.



01

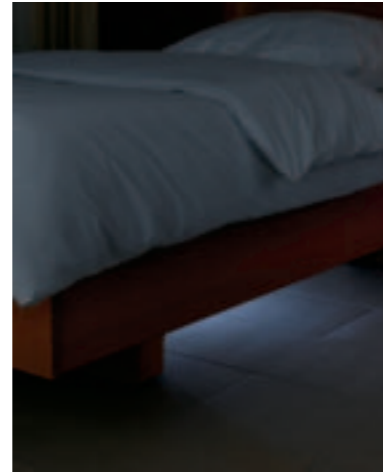

02



03



“No one loves life more than a person who is growing old.”

*Lucius Annaeus Seneca*

Unadorned and contemporary, the contempora M09 care bed most reflects a modern living ambience. The contempora line of bedside cabinets also comes in the attractive wood tone beech 01. The awarded extending bedside table makes everyday life simpler for both residents and care staff. Its one-hand operation allows the bedside table to be easily folded out so that food trays no longer have to be put down.

# contempora M09





# care bed contempora

Safety. Comfort. Functionality.

## 01 Side protections

This innovative safety concept provides major advantages:

- Gapless security
- Can be used as an aid for mobilisation and getting up
- Fall prevention in accordance with the "expert standard Fall Prophylaxis\*" is possible
- Adjusts to four height levels
- Accommodates mattresses up to 18 cm in height for preventing decubitus
- Secure locking mechanisms in accordance with German Federal Institute of Drugs and Medical Devices (BfArM)\*\*
- Alternatively continuous side protections are available

## 02 Hand control

contempora beds are known for their easy operation:

- simple, ergonomic hand control
- Hand control and locking mechanism in one
- Automatic comfort sitting position
- Automated circulation support for residents with vein disease and prevention of decubitus as per Fowler (30° semi Fowler position)

## 03 Comfort bedding area

contempora beds stand for comfort and healthy lying:

- comfort spring technology for better soft-resting effect
- Orientation zone to support sitting and sleeping
- Double retraction of backrest (11 cm) and leg rest (6 cm) for prophylaxis decubitus
- Free choice of mattress
- Removable partitions allow easy cleaning
- Optimal mattress airing



## 04 Easy-click brake

The easy-click system for secure braking:

- Braking and rolling independent of power supply
- Easy-on-the-back operation at any incline height
- No time-consuming adjustment into the lowest position necessary

## 05 Night light

An optional night light can be integrated:

- Easy orientation at night
- More safety for the resident
- Easy power supply control

## 07 Elevation technology

Innovative elevation technology represents improvement and added value:

- 36 -76 cm height adjustment range
- 200 kg safety work load
- Adapted for all lifters
- Free foot space for easier caregiving

## 06 Height adjustment

For care comfort and greater utility:

- Lowest position: 36 cm for easy getting in/out of bed
- Height adjustment: 40 cm in 18 seconds



**Safety:** The Quality and Safety management works in accordance with the regulations DIN EN ISO 9001:2008 and DIN EN ISO 13485:2003 + AC:2007, our beds carry the CE mark. 3 year factory guarantee. **Made in Germany**

**Our recommendation:** "Expert standard fall Prophylaxis" (cf. Expert standard fall Prophylaxis (Prevention) dated Feb. 2005, Chap. 2.6.2.7., p. 71/72, Information can be found at [www.dnqp.de](http://www.dnqp.de))



# secure side railings

**Innovation**  
Version 10


**Comfort Position 1:** with the side protections in their original position, the contempora bed looks very homely.



**Bearing Position 2:** the side protections help residents get their bearings.



**Freedom of move**  
Position 3 at head: this partial position of side safety railings provides safety without overly limiting freedom of move.

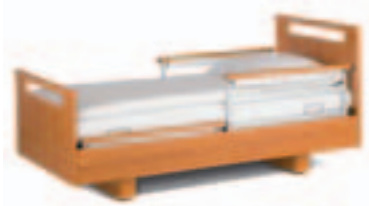




**Security with 12 cm high mattresses**  
Position 3: for safety without any gaps when using 12 cm mattresses.



**Security with 18 cm mattresses**  
Position 4: even with the use of mattresses as high as 18 cm, complete safety is assured.



**Tradition**  
Version 09

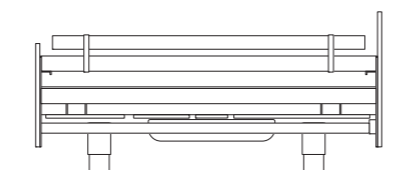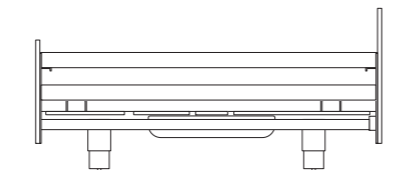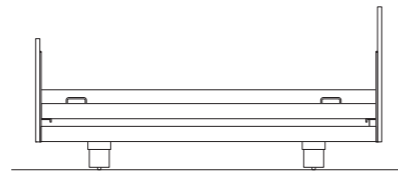


**The security concept** of contempora beds meets the complete range of needs of both residents and nursing staff. The large variety of adjustment options they offer is a major feature, allowing bed security to be adapted flexibly, quickly and as desired to the situation at hand. In their highest position the two sections of the side protections lock together so tightly

that they are every bit as secure as a one-piece full side rail. This makes contempora H10 suitable for any care level including use with restraining belts. In its fourth extension position contempora H10 is also suited for use with special mattresses of up to 18 cm in height. Of course, this care bed is also available upon request with one-piece full side rail.

# mobility



**Left** The bed is movable in any position, allowing great mobility and time savings.  
**Right** contempora with electric braking mechanism is only movable in the lowest position.

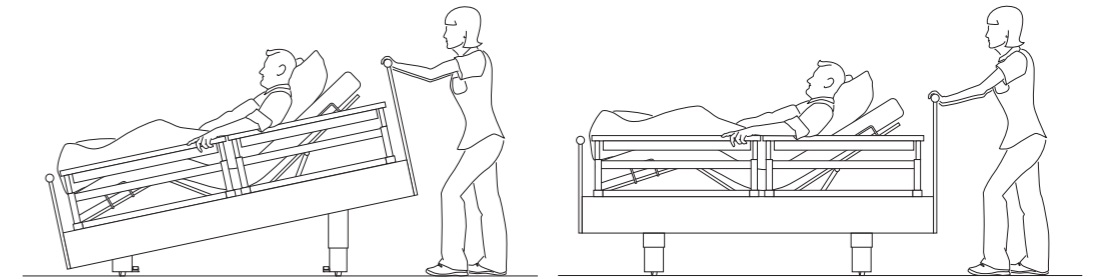


**Easy click brake** The easy-to-reach brake is great for its easy-on-the-back operability and absolute dependability. Easy click provides a large degree of comfort in moving and braking without depending on a power supply.



**Lifting column covers**  
This new elevation technology is hidden inside a slender lifting column cover. The lowest position of the bed at 36 cm makes getting out easier and facilitates transfer into and out of transport chairs. Unoccupied, the height adjustment range of 40 cm is reached within 18 seconds. The bedding area elevates up to a height of 76 cm, making it extremely easy on the back. The stepless height adjustment can be adapted to the body height of individual care staff.

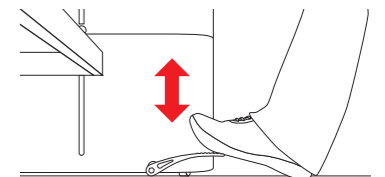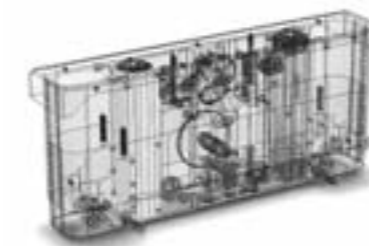

**Secure load capacity**  
The bed guarantees stable security with loads of up to 200 kg.

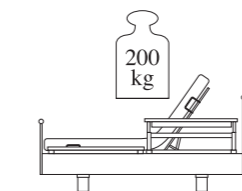

**When you need mobility**, contempora care beds are excellent for their manifold practical advantages. As the beds are movable in any reclining position, no extra trouble is necessary to lower the back support to a flat position. Alongside flexible and spontaneous mobility, the adjustable lifting columns also make work easier: no extra space is required even in pivoted positions. The integrated synchro-function provides for a horizontal reclining surface even with uneven load distribution.




# bedding areas



for triple decubitus propylaxis

## 1

**Comfort spring** The comfort spring technology provides an extra soft resting effect and comfortable reclining – with 12 cm mattresses as well. Harder spring elements along the edge areas create an orientation zone for sitting and sleeping.

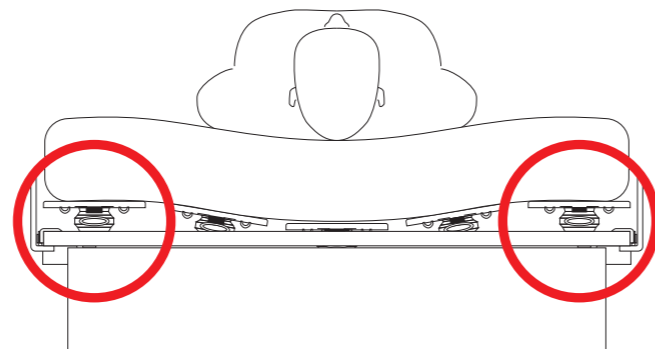


**Comfort** – The comfortable bedding area with its excellent soft resting qualities provides effective prevention of decubitus. The bedding area partitions are removable for easy cleaning.

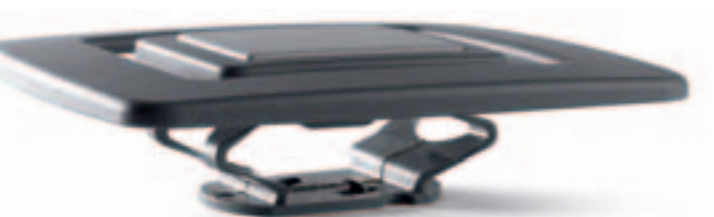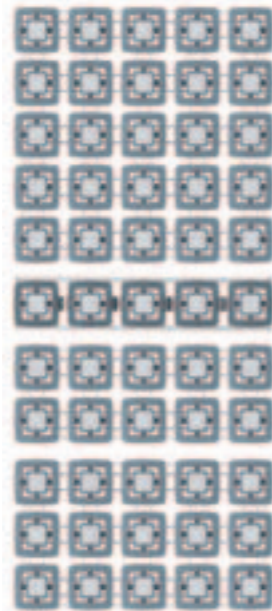

## 2

**Double retraction** Double retraction of back (11 cm) and leg rests (6 cm) decreases pressure on the abdomen area, lowering friction and shearing. Risk of decubitus is reduced.

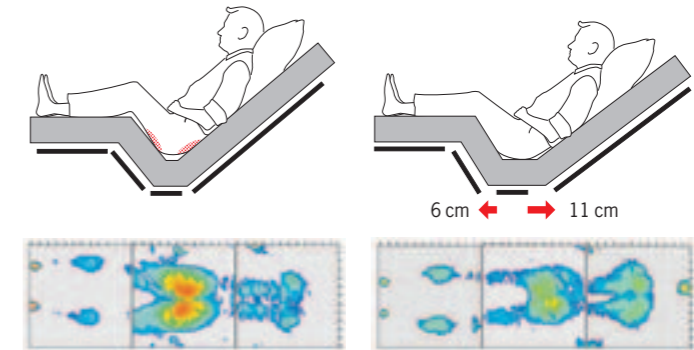


Standard bedding area

**Double retraction** – 17 cm more space in the bed for the auto contour and comfort seating position


## 3

**Free choice of mattress** The system allows effective use of all foam and air cell systems on the market up to 18 cm in height.

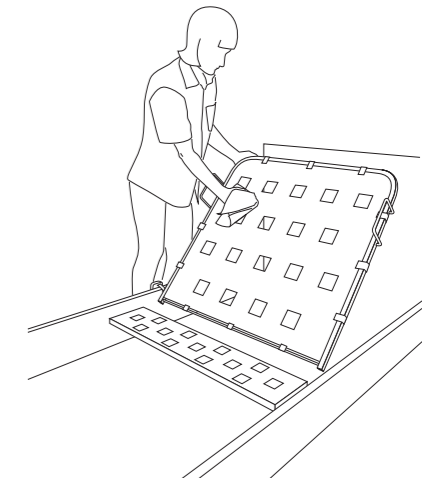

**Left: Aero** – The interlaced wire frame provides for good airing.

**Right: Easyclean** – the ABC covers are especially easy to clean.

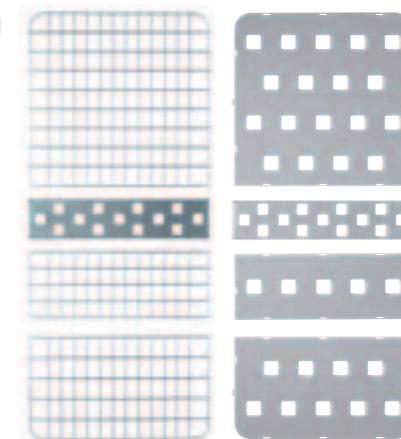

**Many residents** are forced to spend a large part of their time in bed. Having the proper, comfortable position is of great importance to best promote the residents health. contemporary beds offer bedding areas built to provide effective prevention of decubitus. At the same time, nursing staff benefit from the system's easy cleaning.



# storage and mobility

**Easy-care bed adjustment** The automated functions pre-programmed inside the hand control are easy to set.

- Simplicity and safety**  
The hand control's advantages:
- Ergonomic design prevents unintentional button activation
  - 5 selection keys and 2 activation buttons, with LED display of the function activated
  - 3 selection keys are provided for the resident
  - 2-in-1 security package: hand control and locking mechanism in one
  - Non-magnet key for setting in locked, resident or care mode

**Automatic comfort seating**

**Automatic prevention**

**Reset function** After selecting comfort seating or automatic prevention, followed by a downward movement, all lying position adjustments are automatically returned to the zero position.







**Top left** Easy bed entry at 36 cm minimum height setting.

**Top right** Ample space for medically correct 120° positioning, side protections for support and getting bearings.

**Bottom left** contempora supports the effective rope-ladder mobilisation technique.

**Bottom right** Partial side protections allow easy exit from bed.

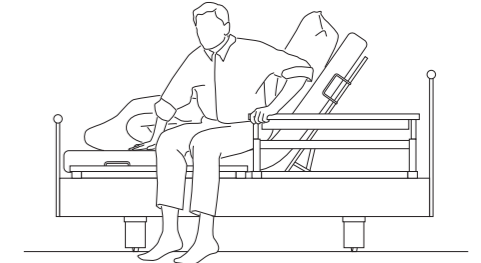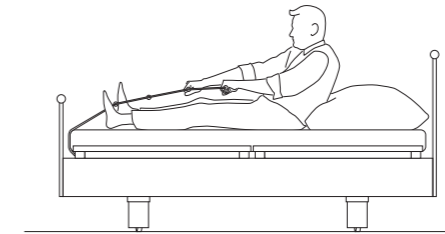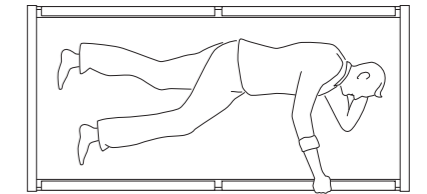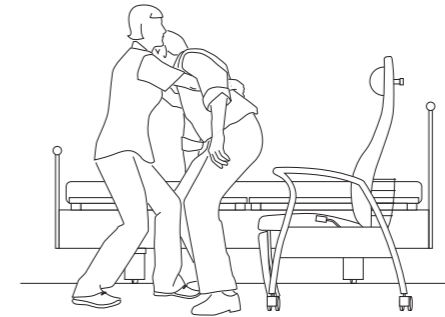


**Automatic comfort seating**

The automatic comfort seating moves the bed to the preset sitting position: convenient for watching television or reading, for eating or communicating at eye level. Mobilisation using the lifter is effectively supported in this position. In addition this position makes breathing easier, as well as digestion. The double retraction in addition provides a relief of pressure in the sacral area.


**Automatic prevention**

supports decongestion with venous diseases and at the same time encourages bed-sore prevention. Through the Fowler-30°-position with the knees raised there is less risk of bedsores. This is achieved through a slightly lowered head position combined with raising the calves, to relieve strain on the heels and joints (e.g. with a cushion placed under the calves). (More information can be found at: [www.decubitus.be](http://www.decubitus.be))

**It is especially when your mobility is decreased** that it is important to be able to rely on dependable technology. That is why contempora beds are equipped so as to provide the highest possible degree of comfort and mobility even for bed-confined residents. Nursing staff also benefit from this, as pre-programmed adjustment functions and caregiver-friendly bed positions make personalised and safe care of residents substantially easier.

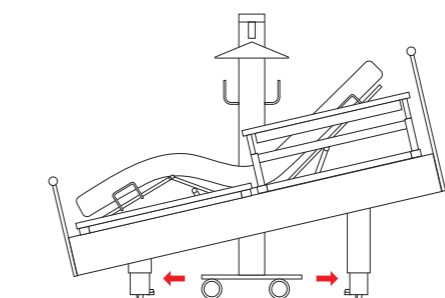

**1. Mobilisation**

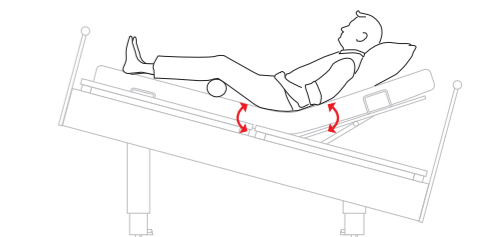

**2. Fowler-30°-Position**  
for bed-sore prevention, based on scientific findings (cf. [www.decubitus.be](http://www.decubitus.be))





# bedside cabinet and server



**The bedside cabinets SD2 and SD3** complete the homely ambience. Their rounded-off corners complement a look that is elegant and modern. The aluminium grips make it easier for residents to open drawers and compartments. A built-in bottle holder also provides for improved handling and organisation.

**Classy design** with rounded-off corners, aluminium grips and glass top plate (optional). **End-to-end aluminium grips** allow effortless pulling out of drawer and compartment. **Compartments** instead of doors: both compartments are easy to reach from the bed and offer lots of storage space.

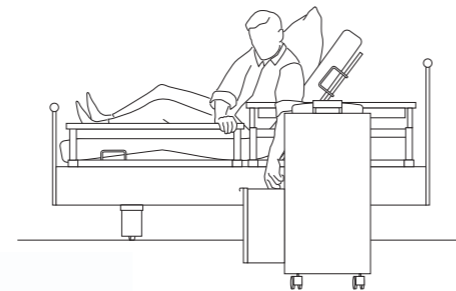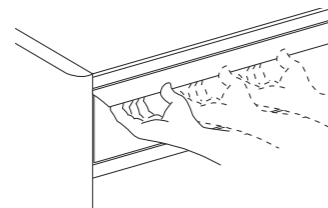


**Combining** the contempora SD2 or vivo E2 bedside cabinet with the stylo S3 server is economical as it allows both products to be used independently – even in separate rooms. The serving stand can be put in park position behind the bedside cabinet SD2 (E2) to save space when not in use.



**The bedside cabinet SD3** is also available with a builtin extending bedside table.

**Convenient** for nursing staff: the extending bedside table, a two-time award-winner for its modern design, is especially easy to position with one hand. The extending bedside table, including height and swivel position, can be adjusted without having to set down the tray. Another plus point: the resident can easily put the bed table in the desired position from bed without assistance.

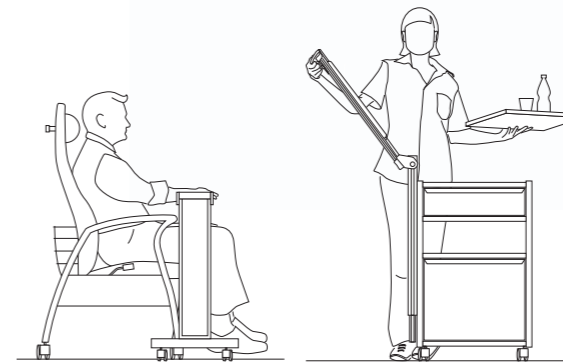


**The U-shaped stylo S3** rolling server stand frame can be used in combination with armrest chairs, transport chairs and wheelchairs. The server stand can also replace as a table.






# product overview

## Side safety protections

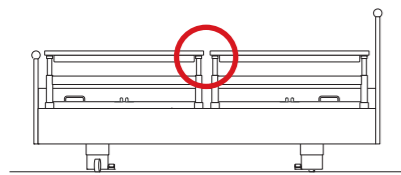


**10** divided side guard, seamless\*  
L 208 cm, W 99 cm, H 37 cm (up to patient surface)

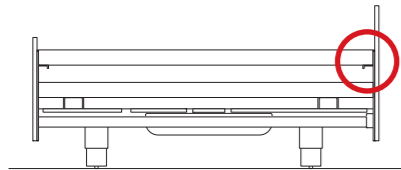


**09** one-piece full side rail  
L 208 cm, W 99 cm, H 37 cm (up to patient surface)

## Head and foot boards


Q

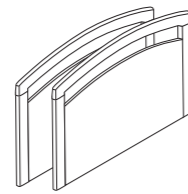


H

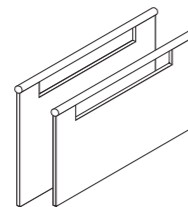


M

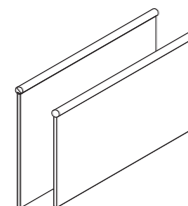

## Moving and braking

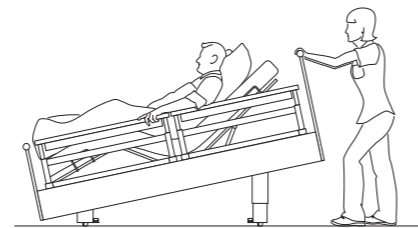


**10-0338** Can be moved and brake set in any incline position

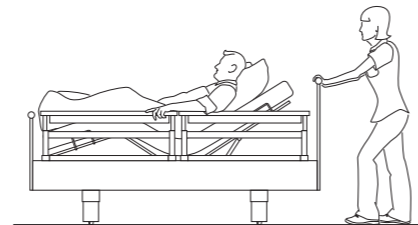


**10-0337** Can be moved and brake set in lowest position

## Night light (Option)

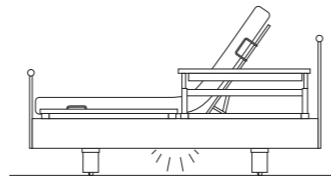


**LED light** under bed helps get night-time bearings


## Bedding area systems

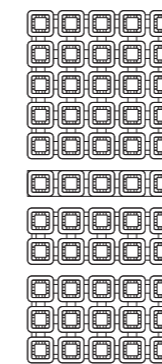

**Komfort**

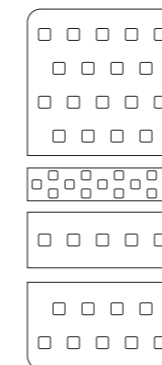


**Easy**

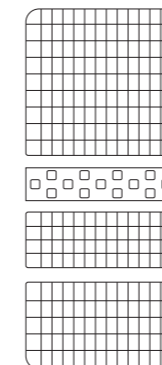


**Aero**


## Bedside cabinet and server

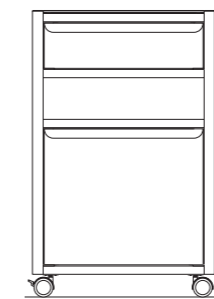

**10-0339** Bedside cabinet SD2  
Glass top plate (optional)

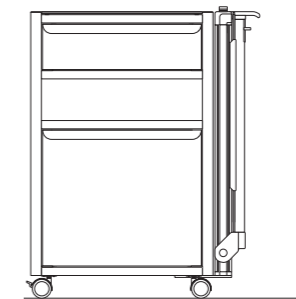


**10-0340** Bedside cabinet SD3  
with integrated extending bedside table glass top plate (optional)

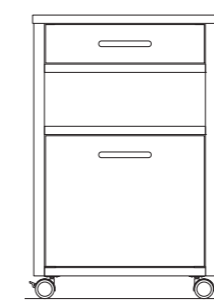

**10-0324** Bedside cabinet vivo E2




**10-0283** Bedside cabinet vivo E3  
with integrated extending bedside table

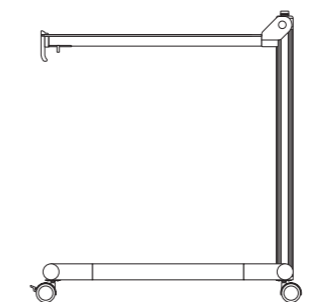


**10-0273** Server stylo S3  
for combining with bedside cabinet without integrated bedside table system


## Options

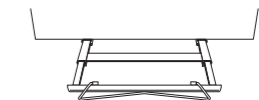

**Linen holder**

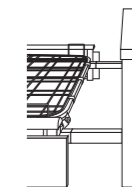

**Bed extension type A**  
side guard 10

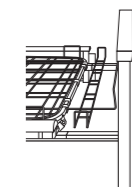

**Bed extension type B**  
side guard 10  
(also for side guard 09, not illustrated)



**Horizontal wall stopper castors**

\* Safe use of patient restraint systems possible in accordance with the BfArM (Federal Institute for Drugs and Medical Devices)




# decors

Concept and Design  
sieger design

Print  
msp druck und medien

Edition  
04/2010, Rev.6, Art.-Nr.: 97-0241  
wissner-bosserhoff

Minor colour variations possible.  
Technical adaptations reserved.

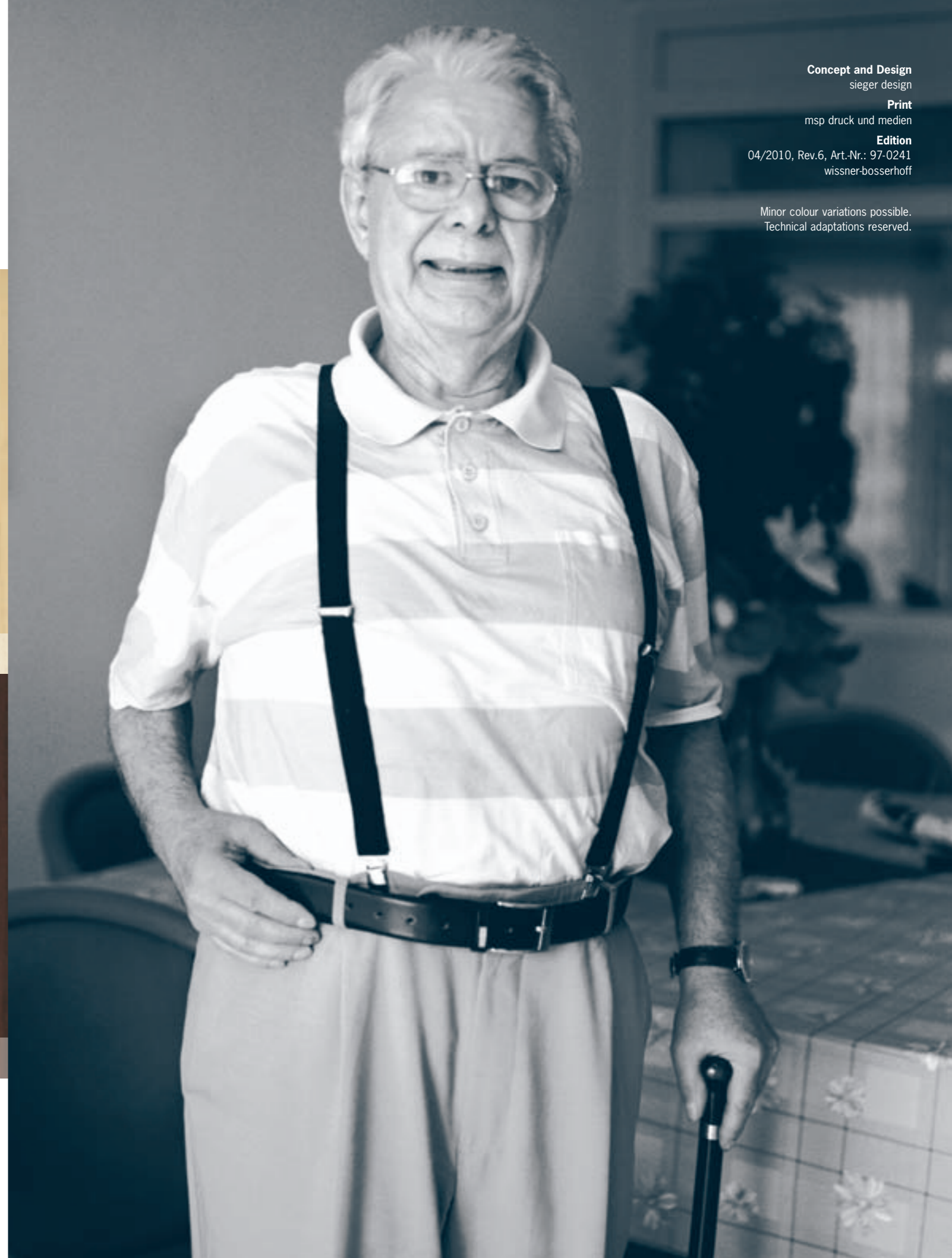

Some illustrated decors available at extra charge.



# wissner- bosserhoff

Furniture for an ambience worth living in

wissner-bosserhoff GmbH  
Hauptstraße 4-6  
D-58739 Wickede (Ruhr)  
Phone: +49 (0) 23 77 / 784-0  
Fax: +49 (0) 23 77 / 784-163  
E-Mail: [info@wi-bo.de](mailto:info@wi-bo.de)

[www.wi-bo.de](http://www.wi-bo.de)

Made in Germany

