

  wissner-  
  bosserhoff

Functional elegance in nursing care

sentida 1-p | 4-p

Psychiatric beds




# Maximum safety and maximum comfort

## sentida – psychiatric beds

Mental and psychosomatic disorders are incredibly costly for the health care system - and this trend is increasing. Mental disorders can seldom be traced back to a single direct cause. Apart from biological factors (e.g., genetics, changes in metabolism) and family predispositions, changing structures in the working world, social pressure, environmental factors, and also democratic change are prime causes for the worldwide increase of depression, schizophrenia, psychosis, and dementia among people from all age groups and walks of life.





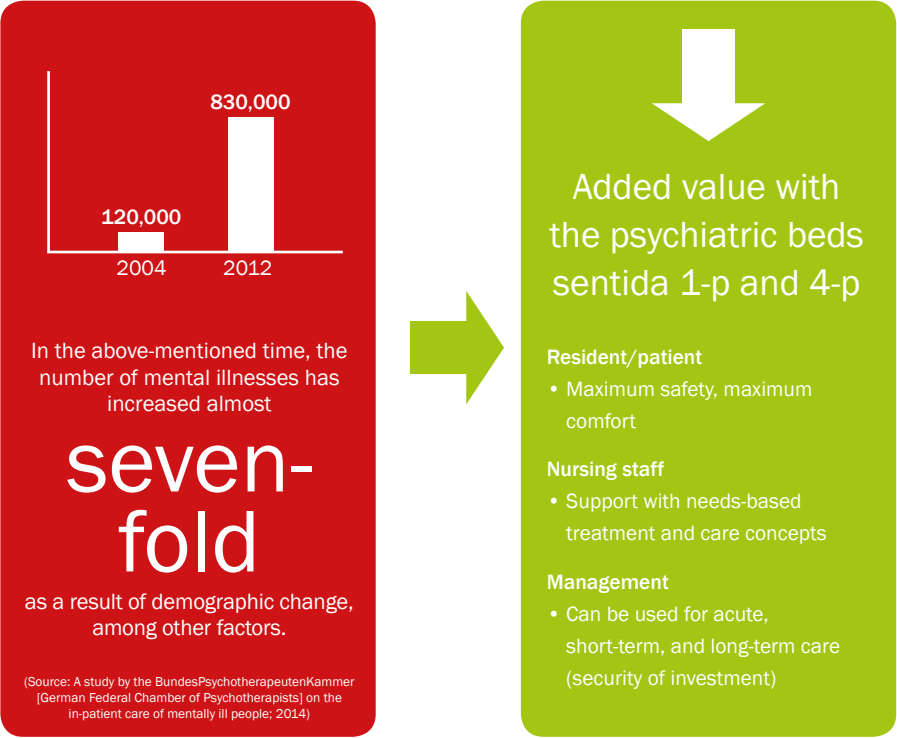



In-patient and out-patient care and treatment concepts in specialized units and institutions for the mentally ill place particular demands on medical devices. sentida 1-p and sentida 4-p are two psychiatric beds that live up to these requirements perfectly. Functional components, for example a non-removable patient surface, integrated restraining belt holders, and built-in, lockable compartments for the handset and power cable reduce the risk of self-mutilation and/or attempted suicide significantly. In addition to the necessary safety aspects in situations of extreme mental crisis, the attractive, homely design of the beds helps the resident/patient get used to the new surroundings.

Specialized for the care of mental patients, sentida 1-p and sentida 4-p are also equipped with the typical functions of a low nursing bed<sup>1)</sup>. These include, among others, the divided SafeFree® side guards and the low position, which contribute to preventing falls when the resident/patient is resting or sleeping.



# sentida 1-p | 4-p

## Portfolio

	 <div>sentida 1-p for perfect stability Item 10-0870</div>	 <div>sentida 4-p for more mobility Item 10-0875</div>
Application	<ul style="list-style-type: none"><li>• Geriatrics, geropsychiatry</li></ul>	<ul style="list-style-type: none"><li>• Clinical psychiatry</li></ul>
Movability	<ul style="list-style-type: none"><li>• In low position</li></ul>	<ul style="list-style-type: none"><li>• At every bed height</li></ul>
Castors	<ul style="list-style-type: none"><li>• Extra-wide feet at both ends</li></ul>	<ul style="list-style-type: none"><li>• 50 or 75 mm double castors at the head and foot end</li></ul>
Braking capacity	<ul style="list-style-type: none"><li>• Electric brake operated using the handset</li></ul>	<ul style="list-style-type: none"><li>• Brake bar at both ends</li></ul>
Safe working load	<ul style="list-style-type: none"><li>• 225 kg</li></ul>	<ul style="list-style-type: none"><li>• 225 kg</li></ul>
Height adjustment	<ul style="list-style-type: none"><li>• 25 – 80 cm (optional: without automatic intermediate stop)</li></ul>	<ul style="list-style-type: none"><li>• 25 – 80 cm (optional: without automatic intermediate stop)</li></ul>
Dimensions	<ul style="list-style-type: none"><li>• Outer dimensions 208 x 106 cm</li></ul>	<ul style="list-style-type: none"><li>• Outer dimensions 208 x 106 cm</li></ul>
Different head/foot sections	<ul style="list-style-type: none"><li>• Low head section + low foot section Accessories: EasyMove 10-0534</li><li>• High head section + low foot section</li></ul>	<ul style="list-style-type: none"><li>• High head section + low foot section</li><li>• Low head section + low foot section</li></ul>
Material	<ul style="list-style-type: none"><li>• HPL head and foot section</li><li>• Side panels optionally either out of HPL or melamine</li></ul>	<ul style="list-style-type: none"><li>• HPL head and foot section</li><li>• HPL side panels</li></ul>
Operating concept	<ul style="list-style-type: none"><li>• Seven-button handset for care with all functions including castor brake release</li><li>• Five-button handset for the resident/patient to adjust the back and leg rests</li></ul>	<ul style="list-style-type: none"><li>• Seven-button handset for care with all functions</li><li>• Five-button handset for the resident/patient to adjust the back and leg rests</li></ul>
Décor	<div><div><div>Volcano white (U11502)</div></div><div><div>Lund birch (R35018)</div></div><div><div>Cottage pine (R55023)</div></div><div><div>Havana cherry (R42006)</div></div><div><div>Natural beech (R24015)</div></div><ul style="list-style-type: none"><li>• Other décors on request</li></ul></div>	<div><div>Volcano white (U11502)</div></div> <ul style="list-style-type: none"><li>• Other décors on request</li></ul>

<sup>1)</sup>For more information on our low nursing beds, please see our brochure "sentida1-6".

# sentida 1-p | 4-p

For special needs in the mental health unit



01 Fixed metal patient surface



02 Individually lockable SafeFree® side guards




03 Individual handset concept

- Lockable handset compartment for nurse handset (all functions)
- Optional resident/patient plug & play handset (not depicted) with limited functions below the patient surface




04 Does not require an external power supply thanks to the optional battery



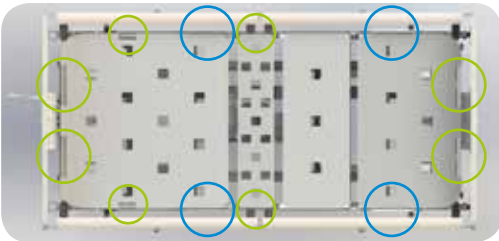



05 Integrated restraining belt holders (○) and additional recesses for additional holders (○)



06 Lockable power cable box in the head section



07 Rubber dampers under the backrest (optional)



08 Cable ducts on the inside; Potential equalization at the patient surface (optional)






# sentida 1-p | 4-p

## Product advantages




### Safety

- The SafeFree® side guards and the handset concept can be tailored to the individual customer needs in in-patient acute, short-term, and long-term care.
- The fixed metal patient surface, covered cable ducts, and integrated restraining belt holders make the bed suitable for exceptional cases in mental health care and increase safety in crisis interventions.

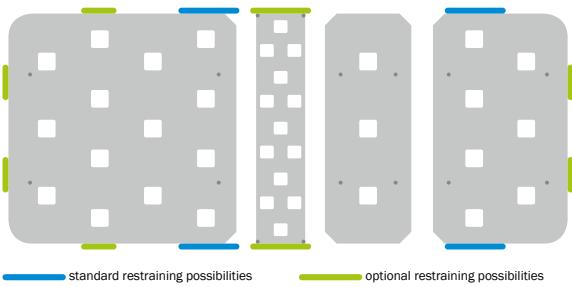

### Improve your processes


- The bed's functionality, design and materials make it possible to use it for the acute, short-term, and long-term care of mentally ill people in nursing homes and hospitals.

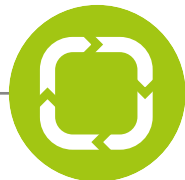

### Improve your workflow

- Three-height concept to adjust the bed to the low, mobilization, and nursing position – promotes resident/patient autonomy and empowerment processes.


**80 cm**  
Optimal nursing position

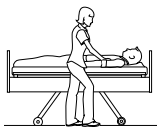

**40 cm**  
Safe mobilization at chair height



**25 cm**  
Fall prevention during sleeping phases

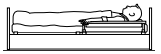

### Hygiene

- Easy to clean surfaces help implement hygiene concepts.
- The scratch and impact-resistant surfaces can withstand manifold strains in daily nursing care.



wissner-bosserhoff voluntarily has its production processes and management audited by TÜV regularly. The company is certified in accordance with the regulations of the standard DIN EN ISO 9001:2008 for quality management systems as well as DIN EN ISO 13485:2003+AC:2007 for the manufacture of medical devices.

---

## Members of LINET Group

### wissner-bosserhoff GmbH

Hauptstraße 4-6 | 58739 Wickede (Ruhr) | Germany

Phone: +49(0)2377.784-0 | Fax: -167 | [info@wi-bo.de](mailto:info@wi-bo.de) | [www.wi-bo.com](http://www.wi-bo.com)



wissner-bosserhoff Germany